

UV-VIS ACCESSORIES

PIKE is addressing the growing need for more sophisticated UV-Vis accessories by featuring research-style specular reflectance accessories, polarizers and automated sampling stages. We also include a selection of most commonly used UV-Vis cuvettes and cuvette holders. Please contact us to discuss customized options.

UV-Vis Cuvettes, Cells, Vials, Holders

Peltier-Controlled Cuvette Holders
*For Experiments Under Tightly Controlled
Temperature Conditions*

Falcon UV-Vis
Precise Cell Temperature Control Accessory

UV-Vis DiffusIR™
Diffuse Reflectance Accessory

UV-Vis Polarizers
Manual and Automated

UV-Vis Spec
Slide Mounted Specular Reflectance Accessories

UV-Vis 10Spec and 85Spec
Fixed Angle Specular Reflectance Accessories

UV-Vis VeeMAX™
Variable Angle Specular Reflectance Accessory

Automated R-Theta Stages
For UV-Vis Spectrophotometers

UV-Vis Cuvettes, Cells, Vials and Holders

PIKE Technologies offers a selection of the most popular cuvettes, vials and cell holders used in the field of UV-Vis spectrophotometry. The cuvettes are manufactured using a heat fusing method which ensures that they are fused into a single homogeneous unit. Cuvettes are carefully annealed to remove any residual strain for maximum cell integrity and physical strength.

The cells can be used with most solvents and acidic solutions. Hydrofluoric acid (HF) and strong bases (pH >9) will negatively affect the cell surfaces.

SPECIFICATIONS

Window Thickness	1.25 mm
Beam Height (Z-dimension)	8.5 mm and 11 mm
Optical Glass Spectral Range	334 to 2500 nm
Far UV Quartz Spectral Range	170 to 2700 nm
NIR Quartz Spectral Range	220 to 3800 nm

ORDERING INFORMATION

STANDARD RECTANGULAR CUVETTES

PATHLENGTH mm	VOLUME mL	OPTICAL GLASS (334–2500 nm)		FAR UV QUARTZ (170–2700 nm)		NIR QUARTZ (220–3800 nm)	
		With PTFE Cover Part Number	With PTFE Stopper Part Number	With PTFE Cover Part Number	With PTFE Stopper Part Number	With PTFE Cover Part Number	With PTFE Stopper Part Number
1	0.40	162-0220	162-0228	162-0236	162-0244	162-0252	162-0260
2	0.70	162-0221	162-0229	162-0237	162-0245	162-0253	162-0261
5	1.70	162-0222	162-0230	162-0238	162-0246	162-0254	162-0262
10	3.50	162-0223	162-0231	162-0239	162-0247	162-0255	162-0263
20	7.00	162-0224	162-0232	162-0240	162-0248	162-0256	162-0264
40	14.00	162-0225	162-0233	162-0241	162-0249	162-0257	162-0265
50	17.50	162-0226	162-0234	162-0242	162-0250	162-0258	162-0266
100	35.00	162-0227	162-0235	162-0243	162-0251	162-0259	162-0267

SEMI-MICRO RECTANGULAR CUVETTES – 4-mm OPENING

PATHLENGTH mm	VOLUME mL	WITH PTFE COVER/CLEAR PART NUMBER	WITH PTFE COVER/BLACK WALLED PART NUMBER	WITH PTFE STOPPER/CLEAR PART NUMBER	WITH PTFE STOPPER/BLACK WALLED PART NUMBER
OPTICAL GLASS (320–2500 nm)					
5	0.70	162-0268	162-0273	162-0278	162-0283
10	1.40	162-0269	162-0274	162-0279	162-0284
20	2.80	162-0270	162-0275	162-0280	162-0285
40	5.60	162-0271	162-0276	162-0281	162-0286
50	7.00	162-0272	162-0277	162-0282	162-0287

FAR UV QUARTZ (170–2700 nm)

5	0.70	162-0288	162-0293	162-0298	162-0303
10	1.40	162-0289	162-0294	162-0299	162-0304
20	2.80	162-0290	162-0295	162-0300	162-0305
40	5.60	162-0291	162-0296	162-0301	162-0306
50	7.00	162-0292	162-0297	162-0302	162-0307

MICRO RECTANGULAR CUVETTES – 2-mm OPENING

PATHLENGTH mm	VOLUME mL	WITH PTFE COVER/CLEAR PART NUMBER	WITH PTFE COVER/BLACK WALLED PART NUMBER	WITH PTFE STOPPER/CLEAR PART NUMBER	WITH PTFE STOPPER/BLACK WALLED PART NUMBER
OPTICAL GLASS (320–2500 nm)					
5	0.35	162-0308	162-0313	162-0318	162-0323
10	0.70	162-0309	162-0314	162-0319	162-0324
20	1.40	162-0310	162-0315	162-0320	162-0325
40	2.80	162-0311	162-0316	162-0321	162-0326
50	3.50	162-0312	162-0317	162-0322	162-0327


FAR UV QUARTZ (170–2700 nm)

5	0.35	162-0328	162-0333	162-0338	162-0343
10	0.70	162-0329	162-0334	162-0339	162-0344
20	1.40	162-0330	162-0335	162-0340	162-0345
40	2.80	162-0331	162-0336	162-0341	162-0346
50	3.50	162-0332	162-0337	162-0342	162-0347

ORDERING INFORMATION

CYLINDRICAL CELLS WITH TEFLON STOPPERS –22-mm DIAMETER

PATHLENGTH mm	VOLUME mL	OPTICAL GLASS (334–2500 nm)	FAR UV QUARTZ (170–2700 nm)	NIR QUARTZ (220–3800 nm)
10	2.80	162-1831	162-1841	162-1801
20	5.60	162-1832	162-1842	162-1802
50	14.10	162-1835	162-1845	162-1805
100	28.20	162-1840	162-1850	162-1810

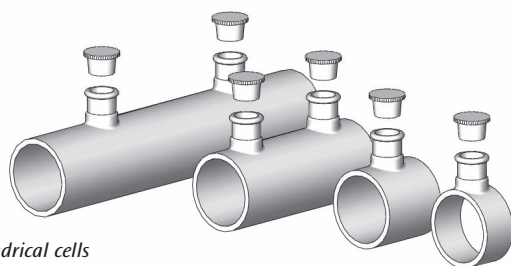

Cylindrical cells

VIALS AND DISPOSABLE CUVETTES


PART NUMBER	DESCRIPTION
162-0205	5-mm Disposable Glass Vials, 5 x 42 mm (200 ea.)
162-0208	8-mm Glass Vials, 8 x 43 mm (200 ea.)
162-0212	12-mm Glass Vials, 12 x 32 mm (200 ea.)
162-0348	Polystyrene, Disposable Cuvettes, 10-mm pathlength, 4.5 mL (100 ea.)
162-0349	Polystyrene, Semi-micro Disposable Cuvettes, 10-mm pathlength, 1.5 mL (100 ea.)

CELL HOLDERS, SPACERS AND STOPPERS


PART NUMBER	DESCRIPTION
111-3650	Cuvette Holder, 10 mm
111-3660	Adjustable Cuvette Holder, 10–100 mm
161-2530	Slide Sample Holder, Cylindrical Cell, 10–20 mm
161-2540	Slide Sample Holder, Cylindrical Cell, 50 mm
161-2550	Slide Sample Holder, Cylindrical Cell, 100 mm
162-1201	Spacer for 1-mm pathlength cuvette
162-1202	Spacer for 2-mm pathlength cuvette
162-1205	Spacer for 5-mm pathlength cuvette

Notes: Please contact PIKE Technologies for replacement Teflon stoppers, covers, and for items not described on this list.


10-mm Cuvette Holder

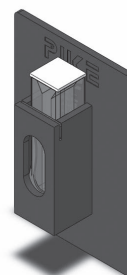


Adjustable Cuvette Holder

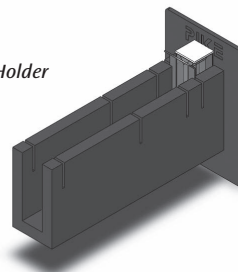

Cylindrical Cell Holder

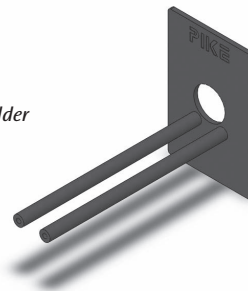

Spacers for 1-, 2-
and 5-mm cells

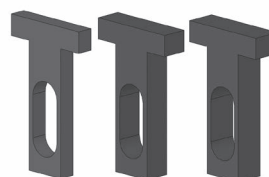

Peltier-Controlled Cuvette Holders for UV-Vis Spectrophotometers – Experiments Under Tightly Controlled Temperature Conditions




FEATURES

- Fast temperature response
- Precision temperature control
- Variable built-in magnetic stirring
- Single- and dual-beam configuration

The new series of PIKE Technologies Peltier accessories includes liquid-cooled cuvette holders, two-channel Temperature Controller and optional Liquid Recirculator. The complete PIKE Peltier Cuvette Accessory requires Peltier cuvette holders, Temperature Controller and Liquid Recirculator. With its efficient design, the Peltier Accessory offers precision temperature control and highly responsive ramping.

The cuvette holder can be used with single- and dual-beam spectrophotometers, and is offered in single or twin configurations. Each temperature-controlled cuvette holder features an efficient heat exchange design, variable speed magnetic stirrers, and a thermoelectric cooler based on Peltier principle. Liquid flow is used to remove excess heat from Peltier elements.


The Peltier cuvette holders are designed to accommodate standard size 10-mm pathlength cuvettes and shorter path cuvettes with appropriate spacers. The holders accommodate 8.5-mm and 15-mm beam Z-height configurations.

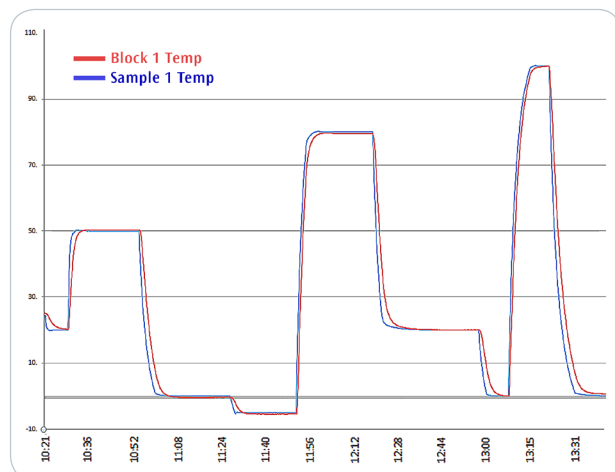


Accessory performance showing block temperature (red) and sample probe temperature (blue).

The temperature controller is configured for a single or dual channel used with the single and twin design, respectively. Within each channel, temperature may be controlled from the precision block RTD sensor or from an optional external temperature RTD probe inserted into the sample cell. All temperature controllers have a USB for PC communication and may be used with PIKE TempPRO software, which provides functions for accessory programming, setting of temperature points and ramping. Temperature controllers have a touchscreen LCD panel

featuring intuitive, menu-driven programming. Single set points or simple ramping functions can be preprogrammed and performed.

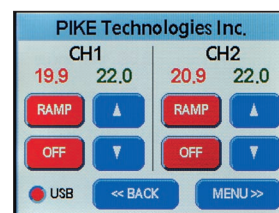

Dedicated temperature controller main menu

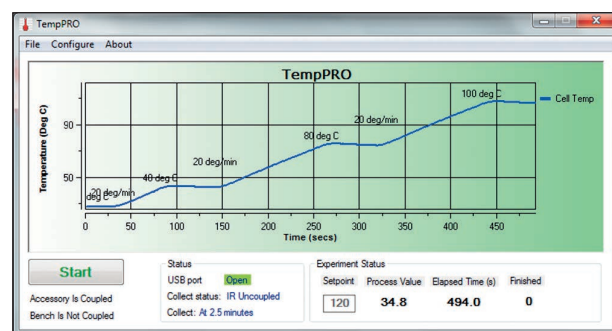

PIKE TempPro software for kinetic experiments

SPECIFICATIONS

Cuvette Holders	
Temperature Range	-5 °C–110 °C
Precision	+/- 0.05 °C
Temperature Accuracy	+/- 0.3 °C from -5 °C to 110 °C
RTD Probe	2 Wire Pt RTD (low drift, high stability)
Stirring Speeds	10 Steps, Variable
Cuvette Size	12.5 mm x 12.5 mm (accommodates smaller size with appropriate spacers)
Z-height	8.5 mm and 15 mm
Temperature Controller	
Precision	+/- 0.1 °C
Channels	Single or Dual
Functions	Set point, temperature ramping
Computer Interface	USB
Input Voltage	90–264 V, auto setting, external power supply
Output Voltage	15 VDC/60 W maximum
Dimensions (W x D x H)	83 x 105 x 85 mm

ORDERING INFORMATION

PART NUMBER DESCRIPTION

171-70XXX	Peltier Cuvette Holder, single – liquid regulated
171-80XXX	Peltier Cuvette Holder, twin – liquid regulated
171-1250	Temperature Control Unit, LCD/PC, single channel
171-2250	Temperature Control Unit, LCD/PC, dual channel
171-1905	RTD Probe for Peltier cuvette
170-1100	Liquid Recirculator for Peltier cuvette
162-1905	Micro Stir Bar, 6.35 mm (l) x 3 mm (dia)

Notes: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)
 Temperature control unit must be selected.

SPACERS

PART NUMBER DESCRIPTION

162-1201	Spacer for 1-mm pathlength cuvette
162-1202	Spacer for 2-mm pathlength cuvette
162-1205	Spacer for 5-mm pathlength cuvette

STANDARD RECTANGULAR CUVETTES

PATHLENGTH mm	VOLUME mL	OPTICAL GLASS (334–2500 nm)		NIR QUARTZ (220–3800 nm)	
		With PTFE Cover PART NUMBER	With PTFE Stopper PART NUMBER	With PTFE Cover PART NUMBER	With PTFE Stopper PART NUMBER
1	0.40	162-0220	162-0228	162-0252	162-0260
2	0.70	162-0221	162-0229	162-0253	162-0261
5	1.70	162-0222	162-0230	162-0254	162-0262
10	3.50	162-0223	162-0231	162-0255	162-0263
20	7.00	162-0224	162-0232	162-0256	162-0264
40	14.00	162-0225	162-0233	162-0257	162-0265
50	17.50	162-0226	162-0234	162-0258	162-0266
100	35.00	162-0227	162-0235	162-0259	162-0267

Falcon UV-Vis – Precise Cell Temperature Control Accessory



FEATURES

- Fast and easy analysis of samples under precise Peltier temperature control
- Choice of cuvette and vial adapters
- Compatible with 1-mm to 10-mm pathlength cuvettes and disposable vials
- Excellent thermal accuracy and precision
- Available for selected UV-Vis spectrophotometers

The PIKE Technologies UV-Vis Cell Temperature Control Accessory is an excellent choice for analysis of liquid samples that require precise temperature control. The accessory is well suited for pharmaceutical, life science and general industrial applications. Temperature range of the accessory is 5 °C to 130 °C with +/- 0.5% accuracy. Heating and cooling is controlled by a built-in Peltier device. The Peltier element provides for reproducible ramping and for reaching target temperatures quickly and reliably. The system is driven by a digital temperature controller – directly, or via PC. Individual sample holders are designed to accommodate standard 1-mm to 10-mm cuvettes (1-, 2- and 5-mm cuvettes require use of spacers) and 5-, 8- and 12-mm glass vials. Sample holders are pin positioned to ensure maximum reproducibility.

The complete Falcon UV-Vis configuration requires the accessory base, vial or cuvette holder, and one of the available temperature controllers. The Falcon accessory is compatible with a selected range of UV-Vis spectrophotometers.



PIKE Liquid
Recirculator

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
110-60XXX	UV-Vis Falcon Base <i>Includes temperature controlled base (digital temperature controller and sample holder need to be selected from the tables below for complete system)</i>

Note: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)

TEMPERATURE CONTROLLERS *(must choose one)*

PART NUMBER	DESCRIPTION
076-1230	Digital Temperature Control Module
076-1430	Digital Temperature Control Module, PC Control (USB)

SAMPLE HOLDERS *(must choose at least one)*

PART NUMBER	DESCRIPTION
111-3610	Vial Holder, 5 mm
111-3620	Vial Holder, 8 mm
111-3630	Vial Holder, 12 mm
111-3640	Cuvette Holder, 1 cm

Notes: Spacers for short pathlength cuvettes are designed to work only with 1-cm cuvette holder.

OPTIONS

PART NUMBER	DESCRIPTION
162-0205	Glass Vials, 5 mm, 5 x 42 mm OD (200 ea.)
162-0208	Glass Vials, 8 mm, 8 x 43 mm OD (200 ea.)
162-0212	Glass Vials, 12 mm, 12 x 32 mm OD (200 ea.)
162-0255	Falcon Quartz Cuvette, 1 cm
162-1201	Spacer for 1-mm Cuvette
162-1202	Spacer for 2-mm Cuvette
162-1205	Spacer for 5-mm Cuvette
170-1100	Liquid Recirculator

SPECIFICATIONS

Temperature Control	Peltier (cooling and heating)
Temperature Range	5 °C to 130 °C
Accuracy	+/- 0.5%
Sensor Type	3 wire Pt RTD (low drift, high stability)
Temperature Controllers	
Digital	+/- 0.5% of set point
Digital PC	+/- 0.5% of set point, graphical setup, up to 20 ramps, USB interface
Input Voltage	90–264 V, auto setting, external power supply
Output Voltage	16 VDC/150 W maximum
Dimensions (W x D x H)	89 x 121 x 83 mm (without FTIR baseplate and mount)

Notes: Peltier device must be water-cooled for proper operation – this is achieved by running cold tap water through the water jacket integrated into the accessory shell, or by the use of an external liquid circulator.

UV-Vis DiffusIR – Highly Efficient Collection Optics for Maximum Sensitivity




FEATURES

- Micrometer-controlled sample focus to optimize results for every sample
- Optional environmental chambers for heating, cooling and high-vacuum applications
- Quick-release feature of environmental chambers for easy insertion and removal of sealed chambers
- Digital PC temperature controller option for programming ramping rates and isothermal hold times

PIKE Technologies introduces the UV-Vis DiffusIR™ diffuse reflectance accessory for research and routine measurements. Powered optical mirrors are diamond-turned aluminum for optimal performance and reflectivity. The UV-Vis DiffusIR base is completely enclosed to shield against external light. To avoid additional stray light from filtering into the accessory, a magnetically-attached light shield is included to cover the sample slide. The standard configuration offers a two-position slide to accommodate reference and sample cups.

The heart of the UV-Vis DiffusIR is a unique monolithic ellipsoidal reflector permanently fixed in place – eliminating the need for repositioning the focus optics for sample placement. The optical design efficiently collects diffuse radiation generated from the sample. The sample Z-position can be optimized by using the micrometer sample focusing adjustment. In this manner the sensitivity of the accessory is maximized without sacrificing precision. The UV-Vis DiffusIR comes equipped with a Sample Preparation and Loading Kit.

Optical geometry of the UV-Vis DiffusIR

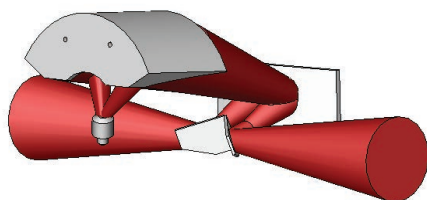


For advanced temperature studies, environmental chambers are available and may be configured for temperatures from -150 °C to 1000 °C. Using the chamber's porous ceramic sample cups, reaction gases may be flowed through the sample. All chambers are compatible with the FT-IR/ NIR DiffusIR also.

Environmental Chamber for the UV-Vis DiffusIR



Coupling the UV-Vis DiffusIR and environmental chambers with the PIKE PC controlled temperature module and TempPRO™ software provides the ability to graphically set up the experiment with up to 20 ramps and hold times.

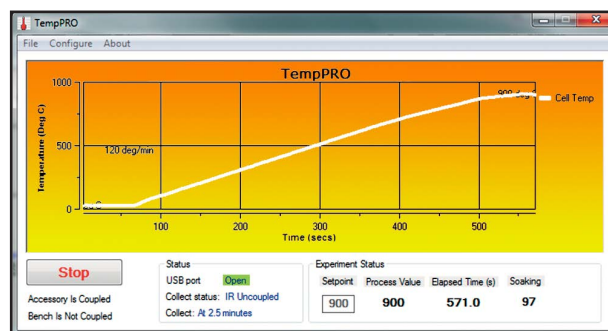

PIKE Technologies TempPRO software provides a graphical interface for temperature control and kinetic measurements.



Liquid nitrogen-cooled system and temperature control module.

ORDERING INFORMATION

PART NUMBER DESCRIPTION

041-10XXX	UV-Vis DiffusIR Accessory <i>Includes Sample Preparation Kit with 2 micro and 2 macro sample cups, sample loading tools, alignment mirror, 35-mm mortar with pestle and KBr powder (100 g)</i>
-----------	---

Notes: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)

OPTIONS

PART NUMBER DESCRIPTION

162-4150	UV-Vis DiffusIR Environmental Chamber, HTV, ambient to 500 °C
162-4200	UV Vis DiffusIR Environmental Chamber, HTV, ambient to 1000 °C
162-4140	UV-Vis DiffusIR Environmental Chamber, LTV, -150 to 500 °C

Notes: HTV and LTV chambers require the selection of a temperature control module. UV-Vis DiffusIR chambers include front plate accommodating environmental chamber (easily changeable with standard UV-Vis DiffusIR front plate), Pin-Loc chamber insertion for easy sample exchange, KBr window, ceramic sampling cups compatible with vacuum and reaction formats, ports and 2 shut-off valves for vacuum operation and ports for connection of water cooling. Operation of the LTV at sub-ambient temperatures requires PN 162-4160 Liquid Nitrogen-Cooled System and Temperature Control Module, and rotary pump for vacuum insulation.

TEMPERATURE CONTROL MODULES

PART NUMBER DESCRIPTION

076-2450	PC Controlled Temperature Module, HTV Chambers <i>Includes digital temperature selection and TempPRO software</i>
076-2250	Digital Temperature Control Module, HTV Chambers
162-4160	Liquid Nitrogen Cooled System and Temperature Control Module, DiffusIR Environmental Chamber, LTV, -150 to 500 °C

Notes: PC controlled temperature module with TempPRO software provides a graphical user interface for setting experiment parameters. Please contact PIKE for PC compatibility. The temperature control modules for the HTV and LTV chambers are not interchangeable.

REPLACEMENT PARTS AND SUPPLIES

PART NUMBER DESCRIPTION

170-1100	Liquid Recirculator
042-2010	Sample Cup, micro, 6-mm diameter, 1.6 mm deep (2 ea.)
042-2020	Sample Cup, macro, 10-mm diameter, 2.3 mm deep (2 ea.)
042-3010	Abrasion Sampling Kit
042-3020	Abrasion Disks, silicon carbide (100 ea.)
042-3025	Abrasion Disks, diamond (50 ea.)
042-3030	Sample Cup Holder and Base
160-8010	KBr Powder, 100 g
042-3040	Sample Preparation Kit
042-3060	Flat Sample Post
042-3080	Alignment Mirror
162-4303	Rotary Pump for vacuum insulation

REPLACEMENT PARTS AND SUPPLIES (cont.)

PART NUMBER DESCRIPTION

160-1132	Disk, KBr, 32 x 3 mm
160-1113	Disk, ZnSe, 32 x 3 mm
160-1372	Disk, UV-Vis SiO ₂ , 32 x 3 mm
160-5049	Disk, SiO ₂ , 32 x 3 mm
160-5125	Disk, low OH SiO ₂ , 32 x 3 mm
160-1143	Disk, CaF ₂ , 32 x 3 mm
162-4210	O-Ring for UV-Vis DiffusIR chamber (10 ea.)
162-4215	O-Ring for UV-Vis DiffusIR chamber cooling line (10 ea.)
162-4251	Ceramic Cup for UV-Vis DiffusIR chamber, porous
162-4270	Alignment Mirror for UV-Vis DiffusIR chamber

Note: Please contact PIKE Technologies for items not described in this list.

SPECIFICATIONS

Optical Design	3X ellipsoidal
Angle of Incidence	30 degrees, nominal
Dimensions (W x D x H)	102 x 225 x 201 mm (excluding light guard tubes and baseplate)
Sample Focus	Micrometer
Sample Positions	2 positions, slide stops for background and sample with no purge loss
Sample Cups	Micro: 6 x 1.6 mm deep Macro: 10 x 2.3 mm deep

ENVIRONMENTAL CHAMBER SPECIFICATIONS

Temperature Range, HTV	Ambient to 500 or 1000 °C
Temperature Range, LTV	-150 to 500 °C
Accuracy	+/- 0.5%
Input Voltage	100–240 VAC (HTV version) 110/220 V switchable (LTV version)
Operating Voltage	28 VDC/84 W (HTV and LTV versions)
Temperature Control	Digital or Digital PC
Heating Rate, Maximum	120 °C/min
Kinetic Setup (requires digital PC controller, includes PIKE TempPRO software)	• Up to 20 temperature ramps • Individual ramp rate and hold time settings • Graphical display of experiment settings • USB interface
Sensor	K Type (for HTV) RTD Type, Pt100 (for LTV)
Vacuum Achievable	1 x 10 ⁻⁶ Torr (13 x 10 ⁻⁴ Pa)
Window Size	32 x 3 mm disk
Leaking Rate	< 6.0 x 10 ⁻¹¹ Pa m ³ /sec
Sample Cup Size	6.0-mm OD, 4.0-mm height 4.7-mm ID, 2.0-mm depth
Sample Cup Design	Porous ceramic compatible with powders and gas flow
Cooling Ports	Quick-Fit, 6-mm ID
Gas/Vacuum Ports	1/8" Swagelok®

UV-Vis Calcite Polarizers – Manual and Automated Versions



FEATURES


- Glan-Taylor and UV Glan-Thompson designs
- High-grade calcite
- High extinction ratio
- Manual versions – 1 degree settable angular resolution
- Automated version – 0.5 degree angular resolution
- Fits in a standard 2 x 3 inch mount

PIKE Technologies offers Glan-Taylor and Glan-Thompson UV-Vis polarizers. These take advantage of the birefringent properties of UV-quality calcite. An air interface is assembled between two right angle calcite prisms in the Glan-Taylor polarizer whereas a UV-transparent cement separates the calcite prisms in the Glan-Thompson polarizer. In both styles, the polarized extraordinary ray passes through both prisms and the ordinary ray is internally reflected and absorbed. The spectral range of these polarizers is 250–2300 nm. Due to the natural origin of calcite the achievable minimum wavelength fluctuates from polarizer to polarizer. However, at 250 nm the transmission throughput is no less than 25%.



Wavelength (nm)	250	300	400	> 500
Minimum Transmission (%)	25	40	65	85

Each polarizer type has a different field of view where the UV-Vis beam is polarized. The UV Glan-Thompson field of view is wider compared to the Glan-Taylor as shown in the figure.


Enhanced angular resolution of 0.5° is gained with the automated version of the PIKE UV-Vis polarizer. This automated feature is advantageous where highly precise angular settings are required and for increasing the measurement simplicity for determining the polarized orientation of a sample. By evaluating the transmission or reflectance of a sample as a function of the automated polarizer angle at a given wavelength, parallel and perpendicular orientation of the sample relative to the polarizer degree setting may be determined. The automated version includes integrated data collection with some commercial UV-Vis spectrophotometer software packages.

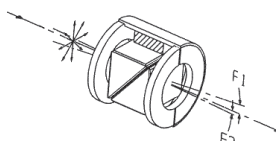

Glan-Taylor field of view

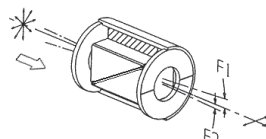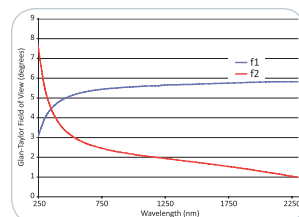

UV Glan-Thompson field of view

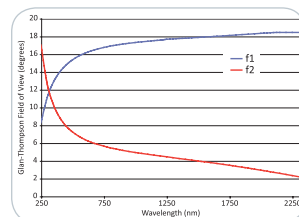

The sample material and sampling configuration dictate the need for a UV-Vis polarizer. Typical materials that may require a polarizer during sampling are often crystals, films, paints, beam splitters, coated glass, and anti-reflective, and anti-glare coatings. Additionally, we recommend using a polarizer for specular reflectance sampling at an angle of incidence greater than 15 degrees where the reflectivity becomes polarization dependent.

SPECIFICATIONS

	Glan-Taylor	UV Glan-Thompson
Material	Calcite	Calcite
Spectral Range	250–2300 nm	250–2300 nm
Clear Aperture	12 mm	14 mm
Extinction Ratio	5×10^{-5}	1×10^{-4}
Manual Polarizer Dimensions (W x D x H)	29 x 50 x 146 mm	49 x 50 x 146 mm
Automated Polarizer Dimensions (W x D x H)	56 x 50 x 146 mm	56 x 50 x 146 mm
Angular Resolution, Manual	1°	1°
Angular Resolution, Auto	0.5°	0.5°

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
198-1623	Manual Glan-Taylor
198-1624	Manual UV Glan-Thompson
198-1625	Automated Glan-Taylor
198-1626	Automated UV Glan-Thompson

Note: Polarizers may not fit in the sample compartments of some smaller spectrophotometers. The automated polarizers include the PIKE Technologies Motion Control Unit and AutoPRO software for automated operation. Please consult PIKE Technologies before placing an order or to inquire about spectrophotometer slide mount holders.

UV-Vis Nanowire Grid Polarizers – Manual and Automated Versions



FEATURES

- Thin profile
- Large acceptance angle, up to 20 degrees
- High transmission and high contrast choices
- Manual and automated versions


PIKE Technologies introduces an innovative new line of high-contrast polarizers covering the UV region and the Vis to NIR region. Using nanofabricating techniques, wire grid lines at a 100-nm pitch are etched on fused silica or glass substrate resulting in a high-performance polarizer. Compared to a traditional calcite polarizer, the large acceptance angle of the nanowire grid polarizer, greater than 20 degrees, eases alignment concerns during use. Additionally, the compact size makes these polarizers ideal for use in confined spaces. The element diameter is 25 mm and has a clear aperture of 19 mm. The polarizer fits a 2 x 3 inch slide mount.

Transmission and contrast ratio of the UV ultra contrast and the broadband polarizers are shown in the next column. Contrast ratio greater than 10,000:1 may be found, making these high-performance polarizers a competitive alternative to calcite polarizers.


There are two manual polarizer types available. The short form has 5 degree scale resolution and the long form has scale resolution of 1 degree. The automated precision polarizers are fully computer controlled and offer the added benefit of increased setting reproducibility with accuracy of +/- 0.5 degree. With automated polarizers an analysis program can be set up through PIKE AutoPRO software, and includes data collection with some spectrophotometer software packages.

SPECIFICATIONS

Substrate	Fused Silica
Performance Range	
UV, Ultra Contrast	240–400 nm
Vis/NIR	300–3200 nm
Element Diameter	25 mm
Clear Aperture	19 mm
Dimensions (W x D x H)	
Manual	50 x 86 x 17 mm
Precision	50 x 146 x 17 mm
Automated	50 x 146 x 55 mm

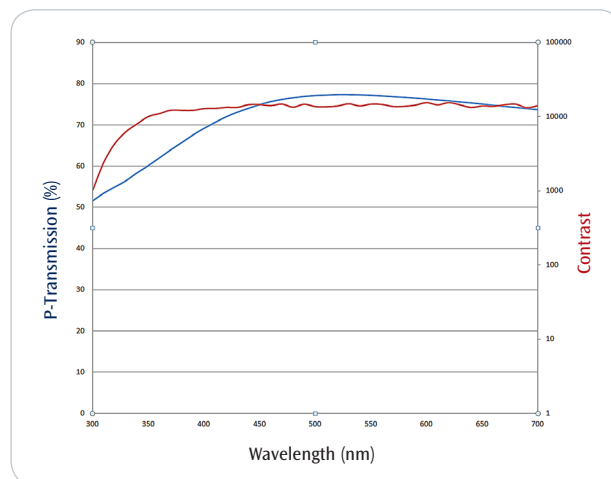

UV Ultra-Contrast Polarizer performance data

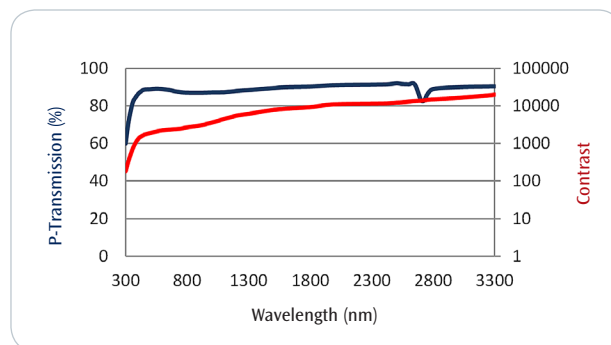


Broadband Polarizer performance data

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
190-2010	Precision UV Polarizer, Ultra Contrast
190-2012	Manual UV Polarizer, Ultra Contrast
190-2015	Precision UV Polarizer, Ultra Contrast, Automated
190-2020	Precision Vis/NIR Broadband Polarizer, High Contrast
190-2022	Manual Vis/NIR Broadband Polarizer, High Contrast
190-2025	Precision Vis/NIR Broadband Polarizer, High Contrast, Automated


Notes: The element diameter is 25 mm. Contact PIKE for 15-mm diameter options. Polarizers may not fit in the sample compartments of some smaller spectrophotometers. The automated polarizers include the PIKE Technologies Motion Control Unit and AutoPRO software for automated operation. Please consult PIKE Technologies before placing an order or to inquire about spectrophotometer slide mount holders.

UV-Vis Spec – Slide Mounted Specular Reflectance Accessories

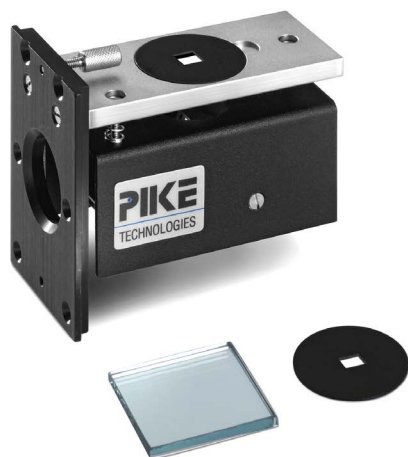

FEATURES


- Fixed angle reflectance accessories including 15, 20, 30, 45 and 60 degrees
- Slide mount design for easy installation
- Optical and anti-reflective coatings testing
- Gloss measurements
- Determining reflectivity of mirrors
- Film thickness calculations

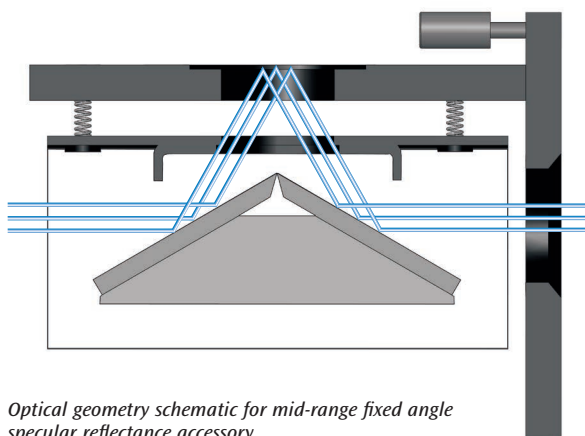

Optical geometry schematic for mid-range fixed angle specular reflectance accessory.

PIKE Technologies relative fixed angle specular reflectance accessories for UV-Vis spectrophotometers span from near-normal to grazing angle. This product sheet highlights our specular reflectance accessories featuring slide mount design. These cover mid-range angles of incidence, between 15–60°, and may be aligned to maximize throughput.

In specular reflectance sampling, the light source is directed to contact the sample at a given angle of incidence. Measured light is collected from the equivalent angle. Typical samples include semiconductors, anti-reflective coatings, color filters, semi-transparent and highly reflective mirrors, optical materials, reflection filters, multiple layers, solar mirrors and solar controlling films on glass. Measured parameters of interest include film thickness, material reflectivity, and coating uniformity and homogeneity.

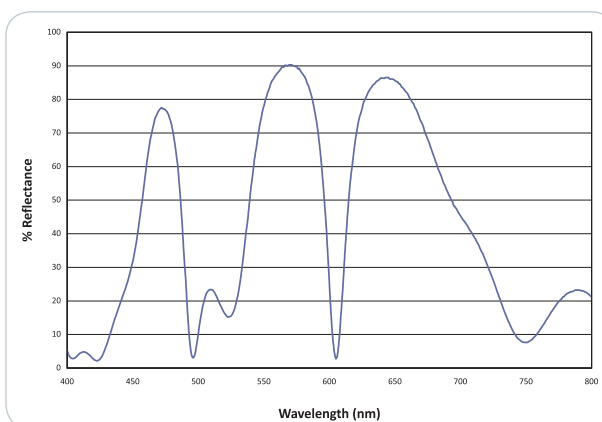

Anti-reflective optics coating spectrum collected using the PIKE UV-Vis 60Spec.

SPECIFICATIONS

Mounting Method	Slide mount
Accessory Dimensions	
Slide Mount (W x D x H)	95 x 51 x 76 mm
Mask Aperture	4 mm x 7 mm
Optical Mirrors	UV-optimized aluminum

ORDERING INFORMATION

FIXED ANGLE SPECULAR REFLECTANCE ACCESSORIES

PART NUMBER	DESCRIPTION
121-1500	UV-Vis 15Spec
121-2000	UV-Vis 20Spec
121-3000	UV-Vis 30Spec
121-4500	UV-Vis 45Spec
121-6000	UV-Vis 60Spec
300-0300	UV Aluminum Mirror (25.4 mm x 25.4 mm)
121-0500	Aperture Mask (4 mm x 7 mm)


Note: Please contact PIKE Technologies to inquire about spectrophotometer sample compartment slide mounts.

UV-Vis 10Spec – Near-normal Sample Reflectivity Measurements



FEATURES


- Measure sample reflectance
- Fixed 10 degree angle of incidence
- Sample illumination using collimated beam precisely fixed at 10 degrees
- Sampling mask sizes of 7, 13 and 25 mm x 4 mm
- Baseplate mount design for stable operation and collection of high-quality spectra

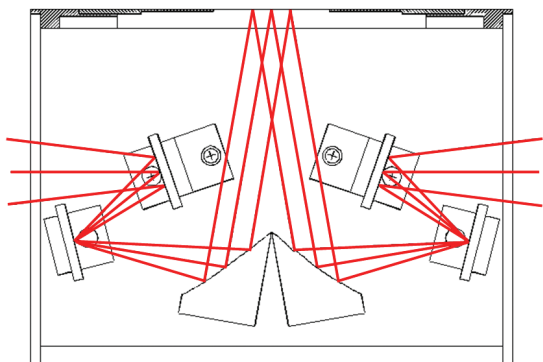

Beam path within the UV-Vis 10Spec specular reflectance accessory.

The PIKE Technologies UV-Vis 10Spec is an optimized specular reflectance accessory designed to make high-performance measurements of sample reflectivity. The UV-Vis 10Spec produces a collimated beam to illuminate the sample area such that the reflectivity measurement is uniformly 10 degrees and not an average of angles produced by a focused beam accessory design. This accessory fits most research-grade spectrophotometers.

The UV-Vis 10Spec may also be used to measure near-normal reflectivity of a wide variety of surfaces including military devices, reflecting optics, anti-reflective (AR) coated surfaces, and other reflecting and non-reflecting materials.

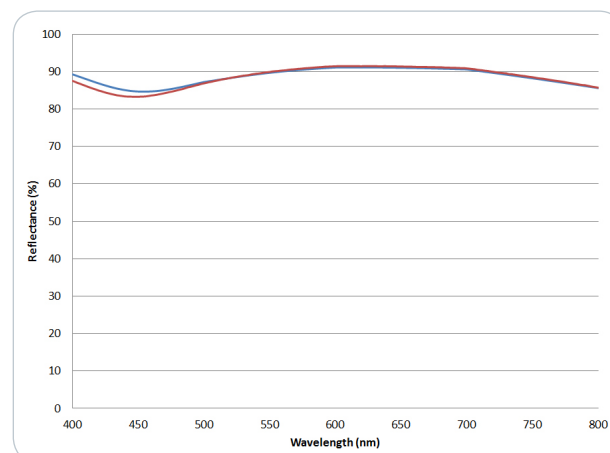

FTIR spectra measuring the reflectivity of SiO_2 coated aluminum mirror with the UV-Vis 10Spec.

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
010-10XXX	UV-Vis 10Spec – 10 Degree Specular Reflectance Accessory Includes 3 sample masks (7, 13 and 25 mm x 4 mm), aluminum alignment mirror and base mount

Note: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)

REPLACEMENT PARTS AND SAMPLING OPTIONS

PART NUMBER	DESCRIPTION
013-4015	Aperture Mask Set, 7, 13, 25 mm x 4 mm
300-0300	UV-Aluminum Mirror, 25.4 mm x 25.4 mm

SPECIFICATIONS


Optics	All reflective
Angle of Incidence	10 degrees
Sample Masks	7, 13 and 25 x 4 mm
Purge Sealing	Purge tubes and purge barb included
Dimensions (W x D x H)	149 x 88 x 118 mm (excludes baseplate)

UV-Vis 85Spec – Specular Reflectance Accessory



FEATURES AND APPLICATIONS


- Small footprint, compact design
- Fixed 85-degree angle of incidence
- Baseplate mount for stable operation and collection of high-quality spectra
- Thin film and coating measurements

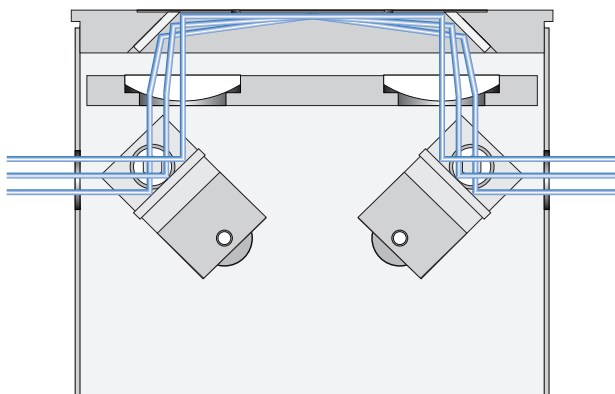

Optical diagram of the UV-Vis 85Spec specular reflectance accessory.

The 85Spec is an 85-degree specular reflectance accessory developed for UV-Vis applications. The accessory is used to analyze mirror-like planar surfaces and thin films on planar surfaces. The grazing angle measurements give the longest possible pathlength through coated samples.

Pre-mounted on a baseplate specific to your UV-Vis spectrophotometer, the 85Spec is an easy-to-use accessory. The mounting allows the 85Spec to be inserted and removed from the sample compartment without alignment.

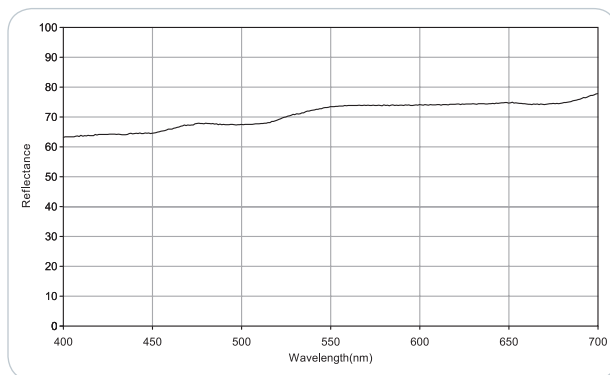

AR-coated ZnSe crystal measured using the 85Spec specular reflectance accessory.

SPECIFICATIONS

Angle of Incidence	85 degrees
Mounting Method	Baseplate
Dimensions (W x D x H)	105 x 84 x 76 mm
Mask Aperture	35 x 3 mm
Optical Mirrors	UV-optimized aluminum

ORDERING INFORMATION

PART NUMBER DESCRIPTION

121-85XXX UV-Vis 85Spec – 85 Degree Specular Reflectance Accessory
Includes an Al substrate alignment mirror and sample mask

Note: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)


REPLACEMENT PARTS

PART NUMBER DESCRIPTION

300-0301 UV-Aluminum Mirror (38.1 mm x 38.1 mm)


121-0501 Aperture Mask (35 mm x 3 mm)

UV-Vis VeeMAX – Variable Angle Specular Reflectance Accessory

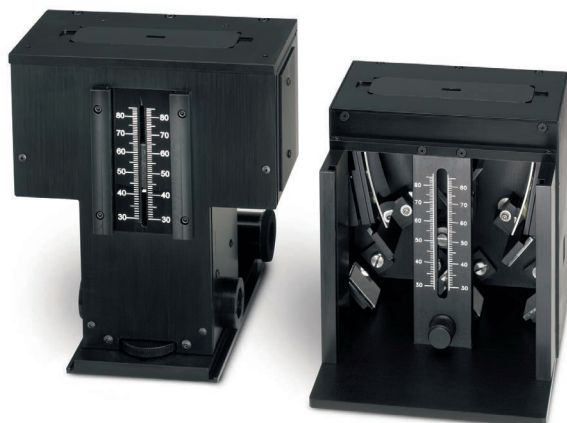

FEATURES AND APPLICATIONS

- Selectable angle of incidence – from 30 to 80 degrees in one degree increments
- Film and coating thickness measurements at optimized angle of incidence
- Grazing angle measurements for the longest path through ultra-thin films
- Characterization of variable reflectivity materials

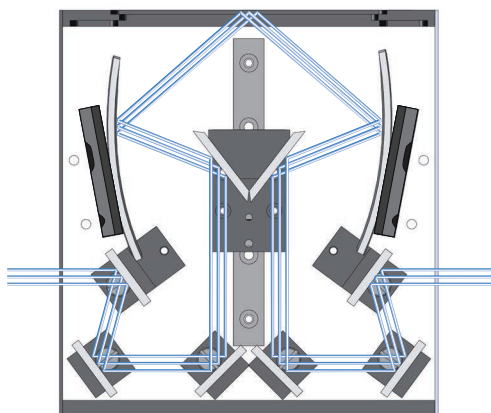


Optical geometry schematic for UV-Vis VeeMAX variable angle specular reflectance accessory.

In specular reflectance sampling, the light source is directed to contact the sample at a given angle of incidence. Measured light is collected from the equivalent angle. Typical samples include semiconductors, anti-reflective coatings, spectral wavelength filters, semi-transparent and highly reflective mirrors, optical materials, reflection filters, multiple layers, solar mirrors, and solar controlling films on glass. Measured parameters of interest include film thickness, material reflectivity, and coating uniformity and homogeneity.

With the UV-Vis VeeMAX™ accessory find versatility for specular reflectance measurements. The angle of incidence from 30 degrees to 80 degrees is easily changed by turning the angle setting dial.

This flexibility allows for optimization of spectral quality for film and coating measurements. In other applications where it is desirable to study the effect of incident angle upon reflected radiation, the UV-Vis VeeMAX offers a sampling solution. For samples such as reflectors, bandpass filters, and hot and cold mirrors. The UV-Vis VeeMAX makes an ideal accessory used to replicate real-world situations such as the effect of a rising and setting sun on absorbance efficiency of solar collectors and architectural glass performance.

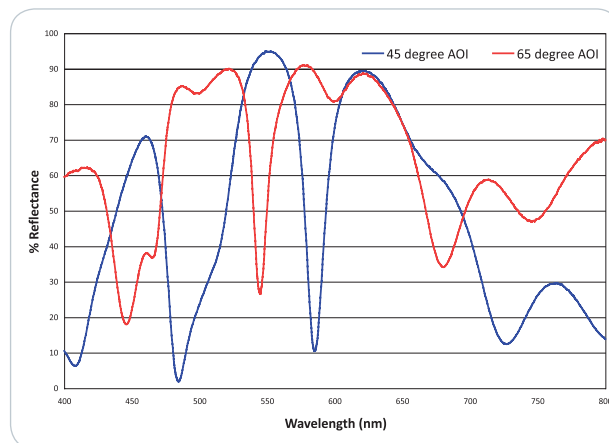

Specular reflectance collected at two different angles of incidence for an anti-reflective coating.

SPECIFICATIONS

Angle of Incidence	30–80 degrees
Mounting Method	Baseplate
Nominal Accessory Dimensions (W x D x H) (dimensions vary with spectrophotometer type)	145 x 135 x 158 mm
Minimum Beam Height	50 mm
Mask Aperture	7 mm x 4 mm 13 mm x 4 mm 25 mm x 4 mm
Optical Mirrors	UV-optimized aluminum

ORDERING INFORMATION

PART NUMBER DESCRIPTION

013-10XXX UV-Vis VeeMAX Variable Angle Specular Reflectance Accessory
Includes masks and aluminum mirror

Note: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)



REPLACEMENT PARTS

PART NUMBER DESCRIPTION

300-0300 UV Aluminum Mirror (25.4 mm x 25.4 mm)

013-4015 Aperture Mask Set (7, 13 and 25 mm x 4 mm)

Out-of-Compartment Microplate Reader – Plate Reading Option for UV-Vis Spectrophotometers

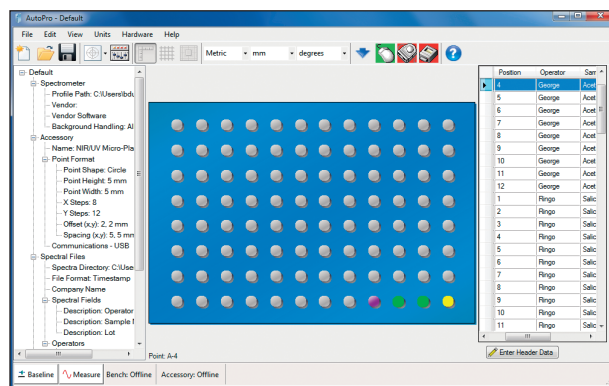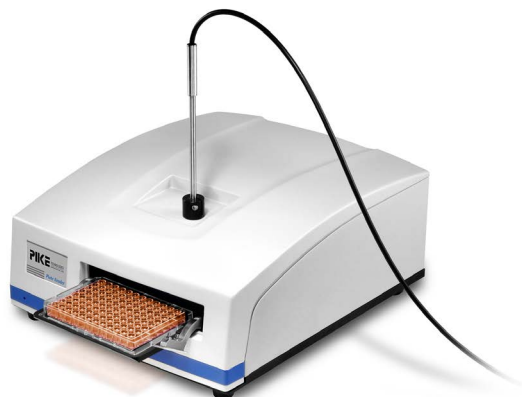

AutoPRO Software – Microplate Reader configuration screen.

FEATURES

- Microplate reading option for standard UV-Vis spectrophotometers
- Scanning and fixed wavelength measurements (spectrophotometer dependent)
- On-board detector
- 6-well to 384-well microplate reading capability
- Custom configurations for automated sampling
- CD-style loading for autoloader interface

The PIKE Technologies Out-of-Compartment Microplate Reader is a unique option available for UV-Vis spectrophotometers. It offers high throughput plate reading capability to a wide range of traditional instruments with standard sample compartments. This allows for flexibility when conducting experiments that require optical configuration for cuvettes, temperature control, integrating spheres (among others) and extra microplate reading functionality when needed. Microplate reading capability is often required in research, drug discovery, bioassay validation, quality control and manufacturing processes in the pharmaceutical and biotechnological industry and academia. The PIKE Microplate Reader can also be adapted to perform automated measurements of filters, optical components and other materials.

The Microplate Reader module features a small footprint and can be positioned next to or above the spectrophotometer. Using a fiber coupler in the sample compartment, the light is sent to the accessory via a single optical fiber probe and collected by an on-board photodiode detector.

The mechanical design of the accessory relies on an X, Y stage with both axes driven by high-precision servo motors with optical encoders for speed and reproducibility. USB and DC power are the only external connections required for this accessory.

Programming and control of the Microplate Reader is done through PIKE Technologies AutoPRO software, which can be integrated easily with many third-party UV-Vis software packages.

SPECIFICATIONS

Optics	Transfer optics module, optical fiber probe and photodiode detector
Accuracy	+/- 25 μ m
Mechanical Specifications	
Repeatability	+/- 5 μ m
Resolution	1 μ m
Run Time	96-well plate – 32 seconds 384-well plate – 84 seconds
Computer Interface	USB
Power Requirements	100–240 Volts AC 50/60 Hz
Dimensions (W x D x H)	11.6 x 13.7 x 6"
Weight	15 lbs

ORDERING INFORMATION

PART NUMBER DESCRIPTION

047-10XXX Microplate Reader Accessory

Note: Replace XXX with your spectrophotometer's Instrument Code. [Click for List >](#)

OPTION

PART NUMBER DESCRIPTION

162-1910 96-Well Polystyrene Sample Plate

REQUIRED OPTICAL COMPONENTS (must select fiber optic probe and photodiode detector)

PART NUMBER DESCRIPTION

047-3011 Photodiode Detector, Cary 50/60

047-3030 Optical Fiber Probe

Note: Contact PIKE Technologies to verify that the accessory is compatible with your spectrophotometer. Requires fiber optic launch optics or fiber couple from spectrophotometer manufacturer.

Automated Transmission R-Theta Rotational Stages for UV-Vis Spectrophotometers




FEATURES

- R-theta motion mapping optical samples
- Complete hardware and software interface package for automated, multi-position measurements
- Light-tight enclosure included
- Custom inserts available

PIKE Technologies' vertical R-theta computer controlled accessories for translating and rotating samples in the spectrophotometer beam. These tools enable transmission mapping of sample surfaces and generating spectroscopy data as a function of sample position. Suitable for determination of film and coating thickness, multilayer film analysis, reflectivity studies and characterization of optical materials. Using the standard sample wheel and custom inserts, these accessories are suitable for analyzing small and large size samples including coated and uncoated glass, optical filters, solar panels and similar materials. Support ring mounts on the accessory's drive and is rotated and translated laterally to produce an R-theta motion covering the entire sampling range of the accessory. Each system incorporates two precision stepper motors for the plate movement.

The operation is managed by PIKE AutoPRO™ software which provides full user programmability and an easy to learn "point-and-click" environment. Polar or X and Y coordinates may be used to define test points. The AutoPRO software allows complex test sequences to be set up, stored as methods and implemented for full flexibility. Each automated accessory is tailored to meet specific sampling needs. This includes adjustments for sample shape/size, and the type of spectroscopic data required. Please contact PIKE Technologies to verify that the selected stage can be integrated into your spectrophotometer's sample compartment and for custom sample inserts.

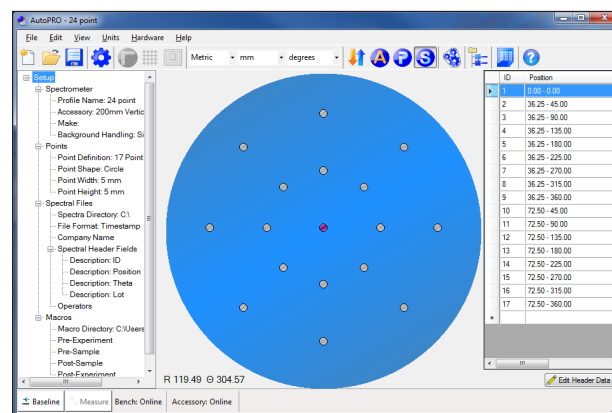

AutoPRO software configured for the Vertical Transmission Accessory.

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
073-4011	Vertical 12 inch Transmission Accessory – Cary 5000/6000 <i>Includes vertical mapping stage with light enclosure, motor controller</i>
073-5100	Vertical 12 inch Transmission Accessory – Lambda 750/850/950/1050 <i>Includes vertical mapping stage with light enclosure, motor controller</i>
073- 5011	Vertical 8 inch Transmission Accessory – Cary 5000/6000 <i>Includes vertical mapping stage with light enclosure, motor controller</i>
073-5010	Vertical 8 inch Transmission Accessory – Cary 500/500 <i>Includes vertical mapping stage with light enclosure, motor controller</i>

Note: For sizes not listed here or custom inserts, please contact PIKE Technologies. Ask about X, Y movement accessories and horizontal stages.

ORDERING TERMS, CONTACT INFORMATION AND GUARANTEE

PART NUMBERS AND PRICE

The PIKE price list includes accessories that may be used with a variety of makes and models of spectrometers. Please specify the part number and description when ordering, including your instrument type and model number. [Click here](#) for a list of spectrometer and spectrophotometer instrument codes. When placing an order, substitute these codes for the final two digits (XX) in the accessory part number.

PIKE Technologies is continually extending the accessory product range. If you are unable to find a required item, please contact us to discuss your needs. We will be glad to assist.

PAYMENT TERMS

Purchase Order Number, cash in advance, MasterCard and Visa are acceptable. Payment is net 30 days, and shipments are FOB Madison, WI USA. Freight charges are prepaid and added to your invoice. If you wish to pay freight charges, please specify this on your order. Prepayment is required for international customers.

INTERNATIONAL HANDLING FEE

For orders placed from outside the United States or Canada, a handling fee of \$40 will apply per order to cover the costs associated with the additional documentation and bank charges required for international shipments.

WAYS TO ORDER

Many products are available for purchase directly through our website. These items are marked on our website with a red shopping cart icon.

Please include the following information when placing an order: your name, phone number, product part number, quantity, ship to address, bill to address, purchase order number and spectrometer model on which the accessory will be used.

Orders may be placed via mail, phone, fax, e-mail or on our website. We accept Visa and Mastercard via phone and direct online purchases. For security purposes, do not send credit card information via e-mail. An electronic order form is available on our website (for P.O. Numbers only – do not use this form for credit card orders). There is no minimum order requirement. Please use the following addresses and phone/fax numbers when placing your orders:

PIKE Technologies, Inc.
6125 Cottonwood Drive
Madison, WI 53719
(608) 274-2721 (TEL)
(608) 274-0103 (FAX)
orders@piketech.com (E-MAIL)
www.piketech.com

DELIVERY

The delivery/shipment date is confirmed upon receipt of an order. Special requirements and custom accessories are subject to different lead times. Please contact us for price quotes and delivery information on these products.

GUARANTEE

All PIKE products are guaranteed to be free from defects in material and workmanship for a period of 12 months from the date of shipment. Should you be dissatisfied, or have any queries, please contact us immediately and we will promptly repair or replace the product at no charge.

PRODUCT RETURNS

Please contact PIKE to receive your Return Material Authorization (RMA) number if you wish to return any of our products. A restocking fee may apply. Customers are responsible for shipping charges for all returned products. For products under warranty, back-to-customer shipping charges will be covered by PIKE. Please do not return any products without obtaining the RMA number first.

TECHNICAL ASSISTANCE

PIKE Technologies offers comprehensive technical assistance. Please contact us via mail, phone, fax or e-mail with your questions.

INTERNATIONAL DISTRIBUTION

PIKE Products are available worldwide. Call or send us an e-mail and we will provide you with an address of the sales office closest to your location. All exports are handled in accordance with the US Export Administration Regulations.

PIKE ON THE WEB

Visit our web site to find out more information about new products, up-to-date PIKE news, pricing, and to see the latest copy of the PIKE Reflections Newsletter! www.piketech.com • info@piketech.com

Customer satisfaction is very important to all of us here at PIKE Technologies, Inc. We have hopefully made the ordering process very fast and easy for you. If you have any questions or concerns about our products or services, please don't hesitate to contact us. We will be happy to make adjustments to fit your needs.

Products and prices are subject to change without notification.

©2017 PIKE Technologies, Inc.

Horizon MB™ belongs to ABB; Luer-Lok™ belongs to Becton Dickinson; Equinox™, IFS™, Quick-Lock™, Tensor™, Vector and Vertex™ belongs to Bruker Optics Inc.; CAB-O-SIL® belongs to Cabot Corporation; Pyrex® belongs to Corning Glass Works; Delrin®, Kalrez®, Teflon®, and Viton® belong to E.I. du Pont de Nemours and Company; Interspec belongs to Interspectrum OU; Spectra Manager™ belongs to Jasco, Inc; Winspec™ belongs to JEOL; EMCOMPRESS® belongs to JRS Pharma; InfraLUM® and SpectraLUM/Pro® belongs to Lumex Ltd; Visual BASIC™ and Windows belongs to Microsoft Corporation; Fluorolube® belongs to OxyChem Corporation; Spectrum™ and Spectrum 10™ belongs to PerkinElmer, Inc; Nujol™ belongs to Schering-Plough; IRPrestige™ belongs to Shimadzu Corporation; Swagelok® belongs to Swagelok Company; Avatar™, Genesys™, Impact™, iS™5, iS™10, iS™50, Magna-IR™, Nexus™, Nicolet™, WinFirst™, OMNIC™ and Protégé™ belong to Thermo Fisher Scientific; Excalibur™, Resolutions Pro™ and Scimitar™ belong to Varian, Inc. All other trademarks are the property of PIKE Technologies.

FTIR AND UV-VIS INSTRUMENT CODES

When ordering a PIKE accessory, replace the **XX** or **XXX** portion of the product's part number with your spectrometer's instrument code below. For assistance, please contact a PIKE customer service representative at (608) 274-2721 or sales@piketech.com.

FTIR INSTRUMENT CODES (XX)

ABB Bomem

FTLA2000-100 (Arid Zone)	80
Michelson 100, MB Series	81
MB 3000	82

Agilent

Excalibur™, Scimitar™, FTS, 600-IR Series	10
Excalibur™, Scimitar™, 600-IR Series with recognition	13

Analect (See Hamilton Sundstrand)

Bio-Rad (See Agilent)

Bruker Optics

IFS™, Vector™, Equinox™ Series.	50
Tensor™, Vertex™ with recognition (Quick-Lock)	51

Buck Scientific

M500	65
------	----

Digilab (See Agilent)

Hamilton Sundstrand AIT

Diamond 20	60
------------	----

Horiba

7000 Series	35
-------------	----

Interspectrum

Interspec 200-X	90
-----------------	----

Jasco

300/600 Series	56
400	57
4000/6000 Series	58

JEOL

Winspec™ Series	46
-----------------	----

Lambda Scientific

Lambda FTIR 7600	66
Lambda FTIR 8600	64

Lumex

INFRALUM FT-02, FT-08	67
-----------------------	----

Mattson (See Thermo Electron)

Midac

M Series	30
----------	----

Nicolet (See Thermo Electron)

Oriel	95
-------	----

Optical Table

	99
--	----

PerkinElmer

1700 Series	70
Spectrum™ GX, 2000	71
Spectrum BX / RX, 1600, Paragon 1000	73
Frontier, Spectrum One, 65, 100, 400 with recognition	74
Spectrum Two with recognition	75

Shimadzu

8300, 8400 Series, IRPrestige™-21, IRAffinity-1s	15
IRPrestige™-21, IRAffinity-1s with recognition (QuickStart)	16
IRTracer™-100	18
IRTracer™-100 with recognition	19

Thermo Electron / Nicolet / Mattson

Infinity, Galaxy, RS Series	20
Genesis™, Satellite, IR 300	21
Impact™ 400, Magna, Protege™, 500 / 700 Series	40
Avatar™, Nexus™, Nicolet™, iS™10, iS™50	40
Model 205/210	41
Nicolet iS™5	42
Avatar, Nexus, Nicolet Series with recognition (Smart)	47

Varian (see Agilent)

UV-VIS INSTRUMENT CODES (XXX)

Agilent/Varian

Cary 50	100
Cary 60	111
Cary 100, 300	110
Cary 4000, 5000, 6000i	120

Jasco

600 Series	600
Optical Table	999

PerkinElmer

Lambda 650, 750, 850, 950 and 1050	700
Lambda 25, 35, 45	730

Shimadzu

1600 and 1700	200
1800 Series	210
2600	240
3600	220

Thermo Fisher Scientific

Evolution 300/600	400
Evolution 200	410

HELPFUL LINKS

[Order Form](#)

[Quote Request Form](#)

[Application Notes](#)

[Visit our Website](#)

[Contact Us](#)



6125 Cottonwood Drive, Madison, WI 53719
(608) 274-2721 (Tel), (608) 274-0103 (Fax)
www.piketech.com · sales@piketech.com